**Jazyky / Languages / Nyelvek**

---

## **THERAPY**

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---

**Operating Instructions (EN)**

 **LIFEFIT®**  
fitness equipment for your life

---

**Dear customers,**

Thank you for purchasing our product. Please read the operating instructions carefully before putting the device into operation.

**1. Technical parameters**

Massage gun with interchangeable heads and the possibility of speed regulation.

Data displayed on the display: speed, battery status

Adapter input voltage: AC 100-240V      Frequency: 50/60Hz

Gun Voltage: DC 24V

Power: 25W

Battery: Li-ion 22.2V, 6x 2000mAh (operating time up to 4,5 hours)

Speed Range (level): 1-20

Amplitude: 9,6 mm

**2. Safety rules**

- The massager is intended for home use.
- If the device does not work properly, stop using it immediately. Contact your supplier.  
Do not use a damaged power cord.
- Do not use the cable in damp areas such as the bathroom, near water, or in water.  
Use in a humid environment may result in electric shock.
- Do not use the device while it is charging.
- Do not touch the space between the head and the body of the gun with your fingers or other objects, as this may cause pinching.
- This product is not a toy. Keep (store) it out of the reach of children.
- People suffering from heart disease, tumors, spinal malformations, people using pacemakers or pregnant women should not use the device.
- The maximum recommended length of one massage is 15 minutes. The device switches off automatically after 10 minutes.

**3. Functions and controls**

The device massages tired muscles, promotes blood circulation or helps with relaxation.

---

(1) "START / STOP"

Turn on the gun with the START / STOP touch switch. Hold the switch until all the numbers are displayed. Use the same switch to turn off the gun.

(2) "SPEED +/-"






Use the touch switches "+" or "-" to increase or decrease the speed to the desired value.

Please touch all touch switches only lightly, there is no need to press them.

**Charging**

The connection for the charging cable is located on the bottom of the gun. Allow the battery to charge for up to 3 hours before first use. When charging, the display lights up and the battery status indicator flashes, when fully charged it lights up continuously. We do not recommend letting the battery drain completely.

**Interchangeable heads**

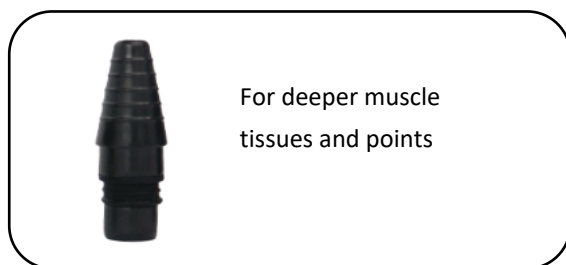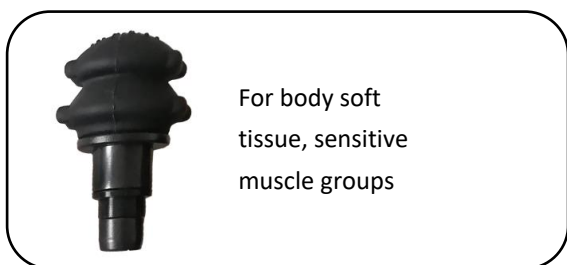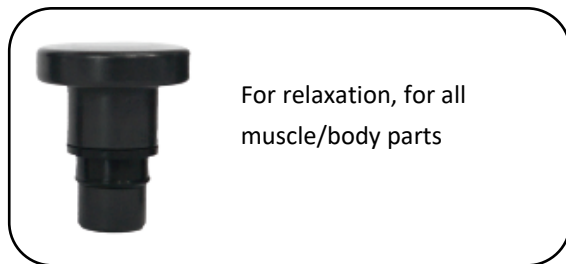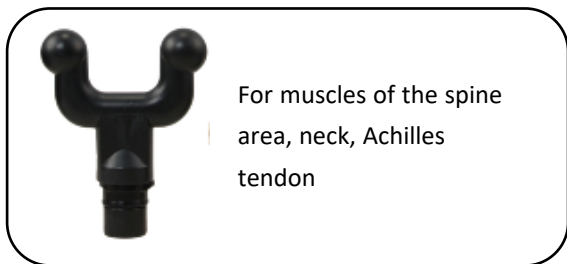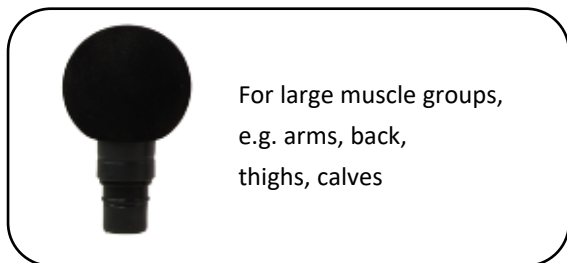---

#### 4. **Maintenance and storage**

- Make sure the gun is unplugged before cleaning. Regularly wipe dust and dirt from the surface of the entire gun with a soft cloth that can be slightly moistened (well wrung out), after cleaning, wipe the gun dry.
- Store the gun in an indoor, dry, clean and well-ventilated area. Protect it from direct sunlight and high temperatures.

#### 5. **Warranty**

The product is protected by the warranty period valid in the country where the sale took place. This warranty covers possible manufacturing defects. It is not provided for wear and tear of the product due to normal use, for defects and wear caused by improper handling, storage and use (e.g. for a purpose other than the product intended) and for any modifications to the product. The sales document serves as a warranty card.

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---

*The author reserves the right not to be responsible for the timeliness, correctness, completeness or quality of the information provided. Claims for liability against the author, which apply to damages of a tangible or intangible nature, caused by the use or non-use of the information provided or the use of incorrect and incomplete information are excluded, if the author does not act in this way intentionally or it is not a gross error.*

*If the purchased (real) product does not completely match the product depicted in the manual, the real product is considered the final and correct product.*

*The photos in this manual may differ slightly in product details from the original.*

Please scan the QR code to view any updated version of the manual.

Verze: 2023/05-1