Kosmetický masážní přístroj na krk a obličej • NÁVOD K OBSLUZE	<b>CZ</b>	<b>3-8</b>
Kozmetický masážny prístroj na krk a tvár • NÁVOD NA OBSLUHU	<b>SK</b>	<b>9-14</b>
Beauty Device for Face and Neck • USER MANUAL	<b>EN</b>	<b>15-20</b>
Arcra és nyakra való masszázrs készülék • KEZELÉSI ÚTMUTATÓ	<b>HU</b>	<b>21-26</b>
Przyrząd do masażu twarzy i szyi • INSTRUKCJA OBSŁUGI	<b>PL</b>	<b>27-32</b>
Gesichts- und Nackenmassagegerät • BEDIENUNGSANLEITUNG	<b>DE</b>	<b>33-39</b>


# Fenité


26/5/2023

*eta*

1

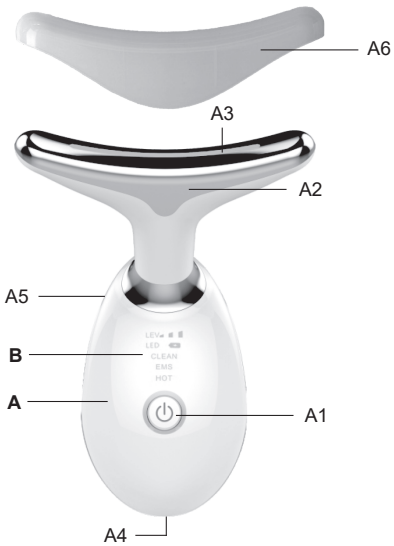B

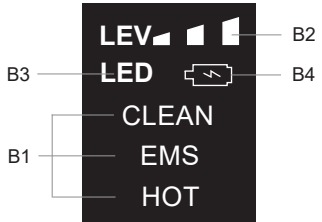C

2

eta 4353

Fenité

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Fenité

## INSTRUCTIONS FOR USE

Dear customer, thank you for purchasing our product. Please keep these instructions including the warranty, the receipt and, if possible, the carton box with the internal packing.

### I. SAFETY WARNING

- Please consider the instructions herein as a part of the appliance and pass them over to the next user of the appliance.
- The appliance is only designated to be used in household and for similar purposes. It is not designed to be used in beauty parlours or for any other commercial utilisation !
- This appliance may only be used by children aged 8 years or older, and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Never immerse the appliance in water!
- When handling the appliance, please proceed so that no personal injury or damage to the appliance will occur.
- Protect the appliance from dust, chemicals, excessive humidity, intensive temperature changes, and do not expose it to direct sunlight.
- Before use, take off all metal items and jewellery (e.g., glasses, watch, earrings, bracelets, rings, etc.).
- Do not insert the appliance or its accessories in any bodily orifices.
- Never use the appliance for other than the intended purposes as described in these instructions for use!
- This appliance contains batteries that are not replaceable.
- This appliance has a heated surface, and the persons with impaired heat perception must be careful in using it.
- **WARNING:** In the case of incorrect use of the appliance that contradicts the operating instructions, there is a risk of injury.
- Texts in a foreign language and images provided on the packaging or on the product, if any, are translated and explained at the end of this language version.
- The manufacturer is not responsible for any loss caused by improper handling of the appliance and its accessories (e.g. injury, harm to health) and is not responsible according to the appliance warranty in the case of non-compliance with the safety instructions above. This appliance is not intended as a substitute for medical care.
- Do not use the appliance if you have hypersensitive, damaged or eczematous skin (e.g. **fresh scars, various wounds, burns, redness, swollen cheeks, skin diseases, facial plastic surgery...**).
- If you perceive pain/ unpleasant feelings or if any other health problems occur, switch off the appliance immediately and consult with your physician or a dermatologist about its use.

- Do not use this appliance if you have metallic implants, metallic dental braces, piercing or an implanted pacemaker.
- Consult the use of this appliance with your physician if you have any of the following contraindications: pregnancy, breastfeeding, allergy, blood coagulation disorder, chronic diseases or if you take sensitive medical preparations sensitive to light.
- The appliance is not intended to cure any skin disease symptoms. Please consult use of the appliance with your doctor, where necessary.
- Do not use the appliance on children, incapacitated persons or persons unable to express their feelings.
- Do not look directly into the light and do not aim at other's peoples eyes.

### NOTICE

#### - The appliance must only be powered by a safe low voltage of 5 V.

- For power supply, use only a USB source complying with valid standards with an output voltage of 5.0 V and an output current of at least 1000 mA (a mobile phone charger with a double-insulated USB connector or a backup power source is suitable).
- The appliance does not have its own power supply.
- Do not use USB ports on other devices (e.g. computer, laptop, tablet, mobile phone, etc.) for power supply. These devices may not be sized for the current consumption of the appliance and could be damaged.

### SPECIAL SAFETY INSTRUCTIONS

- Children must be kept away from the appliance and its adapter.
- Never use the appliance if the adapter is damaged, if the appliance does not operate properly, fell on the ground and suffered any damage, or fell into water. In such case bring the appliance to the specialised service centre to verify the appliance safety and proper operation.
- For additional protection we recommend installing a protective switch (RCD) in the electrical circuitry of the bathroom, with the nominal actuating current not exceeding 30 mA. Ask an inspecting technician or an electrician for advice.
- The appliance and the connected adapter must not be used at places where it could fall into bath tub, wash basin or swimming pool. If the appliance fell into water anyway, do not take it out! First of all, unplug the adapter from the socket and only then take the appliance out.

Then bring the appliance to a specialised service centre to check its safety and proper functionality.

- If the appliance is used in a bathroom, it must be disconnected from the mains after use by unplugging the power cord from the wall socket, because the appliance poses a risk if near to water, even when switched off.
- Always unplug the appliance from charger if you leave it unattended and cleaning.
  - The mains adapter supplied with the appliance is only designated for this appliance.
  - Never plug the adapter in the electrical socket or unplug it with wet hands.
  - Never immerse the adapter into water or other liquids (not even partly) and do not wash it in running water!
  - The appliance gets hot when charging, which is absolutely normal.
  - Charge the battery at a temperature close to the room temperature.
  - Never charge the appliance in direct sunlight or close to heat sources, and if the temperature is below 10 °C or above 40 °C, adversely affects the accumulator service life.
  - Prevent hanging the power adapter cord loosely over desktop edge where children may reach it.
  - The power cord of the adapter must not be damaged with sharp or hot objects or open fire, and must not be immersed into water or bent over sharp edges.

#### Note

ETA a.s. is not a licensed provider of health and cosmetic care and does not have expertise in the field of treatment, diagnostics of the health condition and related effects or suitability for the use of appliances that may affect the health condition in any way. During any exercise, massage therapy or cosmetic intervention, there is a theoretical possibility of physical injury or, for example, worsening of the symptoms of existing diseases. By using this appliance, you agree that you use it at your own risk and assume all risks of injury or deterioration. You also agree to waive any claims (whether known or unknown) against ETA a.s. arising from improper use of this appliance or through negligence. Before using any appliances that may affect your health, first consult your doctor about the suitability of their use.

## II. DESCRIPTION OF THE APPLIANCE (Fig. 1)

### A – Appliance body

- A1 – ON/OFF button + operating mode setting button
- A2 – application head
- A3 – LED (blue/red) for photon therapy
- A4 – USB charging socket
- A5 – touch sensor
- A6 – cover

### B – LED display

- B1 – operation modes (CLEAN / EMS / HOT)
- B2 – LEV symbol of operation intensity (3)
- B3 – LED light indicator
- B4 – Battery level indicator

### C – USB cable

### III. INSTRUCTIONS FOR USE

Before using the appliance, its battery must be charged first. Charging is possible with the enclosed USB cable from the mains, i.e. an AC adapter, or a backup power source.

#### CHARGING

- Before the first use or if the appliance was not used for a long time, charge it completely.
- The standard charging time is **approximately 2 hours** and the fully charged accumulator allows for the operating time of **45-170 min** (based on set regime). The operating time decreases depending on the appliance loading and due to the accumulator ageing and wear and tear.

1. Connect the machine via USB cable **C** to the electrical source (charger). The appliance cannot be used while charging.
2. Do not exceed the designated charging time!

#### Note

The status and charging progress of the accumulator is shown on the display using the **B4** indicator. Blinks while charging and stays on when fully charged. A discharged battery is signaled by the blinking of the **B4** indicator during operation.

#### Optimise the battery service life

- Do not leave USB cable connected to the appliance.
- During charging or regular operation, the appliance heats up, which is a normal phenomenon, and not a reason for a complaint.

#### CONTROL

##### On/off button + setting operation modes (A1)

- By long pushing of button **A1** (approx. for 2 sec.) you switch the appliance on/off. Upon switching on and off an acoustic signal sounds.
- Press short **A1** button to switch over the modes. – The operation mode you have set is shown on display **B** by lighting up its name (CLEAN / EMS / HOT).

#### MODES OF OPERATION AND THEIR USE

The appliance contains 3 different modes with different manifestations. Each works with a different color, intensity of vibration and electricity. pulses (including positive and negative ions), possibly by temperature. To activate the function, it is necessary to hold the body of the product (thereby touching the **A5** belt sensor and closing the loop). We recommend that the process be carried out gradually as shown in Fig. 2.

**Each mode can only be used once a day, at the same time we recommend using all three modes consecutively within one treatment cycle.**

#### Fotons therapy

**Red light** – effectively penetrates deep into the skin, revitalizes collagen, softens wrinkles, reduces pigmentation, increases the efficiency of absorption of cosmetic products, improves blood circulation in the subcutaneous tissue.

**Blue light** – it facilitates the synthesis of collagen, rejuvenates the skin, removes sagging skin, softens skin pores, has an antibacterial and anti-inflammatory effect, and contributes to soothing the skin (most suitable for oily, sensitive and acne-prone skin).

**Green light** – supports standard cell function, accelerates lymph circulation, reduces swelling and soothes the skin.

**Note**

- **Each mode can only be used once a day, at the same time we recommend using all three modes consecutively within one treatment cycle** (individuals with more sensitive skin 2-3 times per week).
- Procedure time approx. 4 minutes.

**CLEAN**

Cleanses impurities in deep pores, shrinks pores, helps reduce pigmentation spots and acne, and contributes to solving the problem of black spots and dull skin color.

**Postup:**

1. Turn on the appliance using the **A1** button, the **CLEAN** mode is the default and is automatically activated. Hold the appliance so that you touch (hand/fingers) the **A5** sensor.
2. Gently massage the desired areas according to Fig. 2.

**EMS**

This mode is intended for easier incorporation of cosmetic products into the deeper layers of the skin. Thanks to this regime, you can effectively apply products such as: moisturizing and rejuvenating creams, hyaluronic acid, oligopeptides, etc. to the skin. It helps to calm the skin.

**Procedure:**

1. Apply the desired cosmetic product to the skin.
2. Turn on the appliance using the **A1** button, set the **EMS** mode and hold the appliance so that you touch (hand/fingers) the sensor **A5**.
3. According to Fig. 2, trace the individual parts of the face/neck. As soon as you feel resistance on the application head, it means that the product has been absorbed into the skin and you can move on to other parts of the face.

**HOT**

This mode is intended for smoothing and softening the contours of the face. Thermal massage supports blood circulation in the skin, increases its elasticity and freshness and keeps the skin in good condition for a long time

**Procedure:**

1. You can apply a facial mask to your face before use.
2. Turn on the appliance using the **A1** button, set the **HOT** mode and hold the appliance so that you touch (hand/fingers) the sensor **A5**.
3. Gently massage the desired areas according to Fig. 2.

**Automatic shutdown function**

During operation, after approx. 4 minutes, an audible signal will sound and the inscription showing the working function (**CLEAN, EMS, HOT**) will start flashing, and after approx. 1 min of operation, the appliance will be switched off.

**IV. MAINTENANCE**

The appliance must be cleaned after every use! Use only soft cloth to wipe the surface of the appliance and the USB cable. Do not use any abrasive or aggressive cleaning agents!

**Guarantee of 24 months does not apply to decrease in the accumulator's effective capacity caused by its use or age. Effective capacity gradually decreases depending on the method of use of the accumulator.**

**Expert service must do the large maintenance or to intervene into inside parts of the appliance!  
Non-observance of the instructions of the manufacturer cancels the right for guarantee repair!**

## V. ENVIRONMENT

If the dimensions allow, there are labels of materials used for production of packing, components and accessories as well as their recycling on all parts. The symbols specified on the product or in the accompanying documentation mean that the used electrical or electronic products must not be disposed of together with municipal waste. For proper disposal, hand them over at special collection places where they will be accepted free of charge. Suitable disposal of the product can help to maintain valuable natural resources and to prevent possible negative impacts on the environment and human health, which could be the possible consequences of improper waste disposal. Ask for more details at the local authorities or at a collection site. Fines can be imposed for improper disposal of this type of waste in agreement with the national regulations. Dispose of the battery safely by means of collection networks intended for this purpose. Never dispose of the batteries by burning!

### Taking out the battery

Only take out the battery from the appliance after it is fully discharged, and after the end of the appliance's service life! In the course of operation, the battery does not need to be removed!

## VI. TECHNICAL DATA

Voltage (V)	Specified on the rating place
Appliance weight (g)	100
Appliance protection class	III.
Dimensions (H x D x W) in mm	50 x 120 x 90

**The manufacturer reserves the right to change the technical specifications and accessories for the respective models.**

**HOUSEHOLD USE ONLY. DO NOT IMMERSE IN WATER OR OTHER LIQUIDS.**

**TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.**

The  symbol indicates a **WARNING**.

**Do not use the appliance close to bath tubs, showers, bath basins or other water containers.**

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© DATE 26/05/2023